


# SELECTION GUIDE

## Dispensers/Dilutors



OUR BEST DISPENSERS!

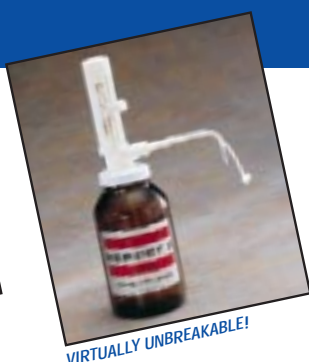

VIRTUALLY UNBREAKABLE!



OUR BEST ALL GLASS DISPENSERS!



OUR LOWEST PRICED DISPENSERS!

### REPIPET DISPENSERS

	Product	Dispenser Volume (ml)	Composition	Reagent Compatibility	Autoclavable	Page
<i>OUR BEST DISPENSERS!</i>	Repipet® Dispensers	0.5, 1, 5, 10 20, 50, 100	Glass Barrel Glass Plunger TEFLON Tubing	All reagents except Hydrofluoric Acid "AL" models not for Organic Solvents	Yes	19
<i>OUR BEST MINI-DISPENSERS!</i>	Mini-Repipet Dispensers	0.5, 1, 2.5, 5, 10,	Glass Barrel Glass Plunger TEFLON Tubing	All reagents except Hydrofluoric Acid	Yes	20
<i>ALL GLASS DISPENSERS!</i>	Glass Joint Repipet Dispensers	0.5, 1, 5, 10, 20 50, 100	Glass Barrel, Glass Plunger, TEFLON Tubing, GR model-Glass Inlet Tubing	All reagents except Hydrofluoric Acid	Yes	21
<i>VIRTUALLY UNBREAKABLE!</i>	Repipet II Dispensers	5, 10, 20	Glass Barrel Glass Plunger TEFLON Tubing HALAR Sleeve	Most reagents except Hydrofluoric Acid & Organic Solvents	No	22
<i>OUR LOWEST PRICED DISPENSERS!</i>	Fixed and Adjustable Volume Dispensers	2, 5, 10	Polyethylene Polypropylene EPR	Non-corrosive reagents	No	23
		0.1, 0.2, 0.5, 0.7, 0.9, 1, 1.5, 2, 2.5, 3, 4, 5, 10	Polyethylene Polypropylene EPR	Non-corrosive reagents	No	23

### REPIPET DILUTORS

		Min. Ratio Diluent to Aspirated Volumes	Composition	Reagent Compatibility	Autoclavable	Page
<i>COMPARE THE PRICE TO AUTOMATIC DILUTORS!</i>	Repipet Single Stroke Dilutors	1:1, 13:1, 2.3:1 4:3	Glass Barrel Glass Plunger TEFLON tubing	All reagents except Hydrofluoric Acid "L" model is not for organic solvents	Yes	24

### PIPETS


		Pipet Volume (ml)	Composition	Page
<i>CHEMICALLY RESISTANT!</i>	Air Displacement Pipets	5, 10	All Plastic	25

# Dispensers/Dilutors

REPIPET® DISPENSERS

Labindustries

Standard and Low Silhouette  
Dispensers

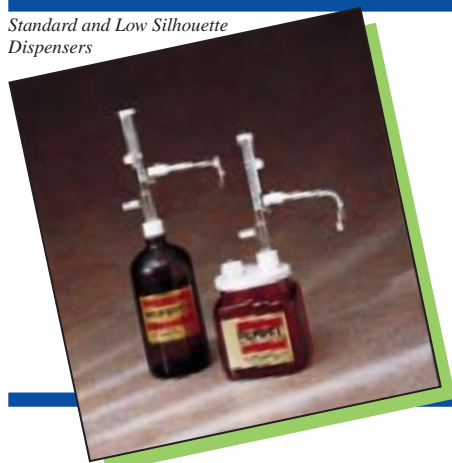

## OUR BEST DISPENSERS!

- — Excellent for use with acids and solvents
- — Easily maintained
- — Reliable performance
- ★ — Glassware carries a 50% lifetime guarantee against breakage

## ORDERING INFORMATION

Model #	Size (ml)	US List Price
<b>Universal REPIPET Dispensers</b>		
3001/2AU	0.5	\$274.00
3001AU	1.0	\$274.00
3005AU	5.0	\$280.00
3010AU	10.0	\$280.00
3020AU	20.0	\$301.00
3050AU	50.0	\$366.00
3100AU	100.0	\$457.00
<b>Low Silhouette REPIPET Dispensers</b>		
3001/2AL	0.5	\$281.00
3001AL	1.0	\$270.00
3005AL	5.0	\$286.00
3010AL	10.0	\$286.00
3020AL	20.0	\$295.00
3050AL	50.0	\$371.00
<b>Standard REPIPET Dispensers</b>		
3001/2A	0.5	\$270.00
3001A	1.0	\$270.00
3005A	5.0	\$274.00
3010A	10.0	\$274.00
3020A	20.0	\$280.00
3050A	50.0	\$348.00

## PRODUCT DESCRIPTION

- Dispenser components are fully autoclavable.
- Full range adjustability, 100-Division Scale.
- 1% accuracy and 0.1% reproducibility at full scale guaranteed.
- Includes a magnifying indicator for easy-reading accuracy.
- Pyrex® and Teflon® fluorocarbon construction for safe use with all liquids except hydrofluoric solutions.
- Models ending in A and AU include both magnifying indicator and metal indicator.

## UNIVERSAL REPIPET DISPENSERS

- Dispense from any screw-cap container with an inner diameter opening greater than 19 mm or any container with a 24/40 or 29/42 standard taper by adding an adapter for a ground-glass joint.
- Supplied with 3-cap and washer sets to fit your container: 28-430, 33-430 & 38-430.
- 45 mm cap & washer set optional. #CV830X1B
- 24" (61 cm) of Teflon tubing to fit your container.

## LOW SILHOUETTE REPIPET DISPENSERS

- Ideal for shelf storage, fits easily in the refrigerator.
- Square bottle reduces the risk of damage from tipping.
- A second opening in the Halar®-fluorocarbon cap allows you to add reagent without disturbing the dispenser.
- Supplied with 1000 ml square amberized bottle.

## STANDARD REPIPET DISPENSERS

- Perfect all purpose choice for dispensing all reagents except hydrofluoric acids, including concentrated acids, concentrated alkalis, chlorinated hydrocarbons, and organic solvents.
- Supplied with a 950 ml round amberized bottle.

48 Hour or Sooner  
Express Shipping  
Guaranteed!\*



\*see inside front cover for details

# Dispensers/Dilutors

## MINI-REPIPET® DISPENSERS

Labindustries

Model MINI3005A Dispenser



### OUR BEST MINI DISPENSERS!

- — For smaller reagent sample volumes
- — Safe for all reagents except HF
- — Air intake filter prevents airborne contamination
- ★ — Glassware carries a 50% lifetime guarantee against breakage

### ORDERING INFORMATION

Model #	Size	US List Price
MINI3001/2A	0.5 ml	\$255.00
MINI3001A	1.0 ml	\$255.00
MINI30021/2A	2.5 ml	\$255.00
MINI3005A	5.0 ml	\$255.00
MINI3010A	10.0 ml	\$255.00

### PRODUCT DESCRIPTION

- Full range adjustability, 100-Division Scale.
- 1% accuracy and 0.1% reproducibility at full scale guaranteed.
- Includes both magnifying indicator and metal indicator.
- 120 ml round, amberized, glass bottle.
- Support stand provides stability.
- Autoclavable.

48 Hour or Sooner  
Express Shipping  
Guaranteed!\*




\*see inside front cover for details

# Dispensers/Dilutors

GLASS JOINT REPIPET® DISPENSERS

Labindustries

Model 3100GR Dispenser

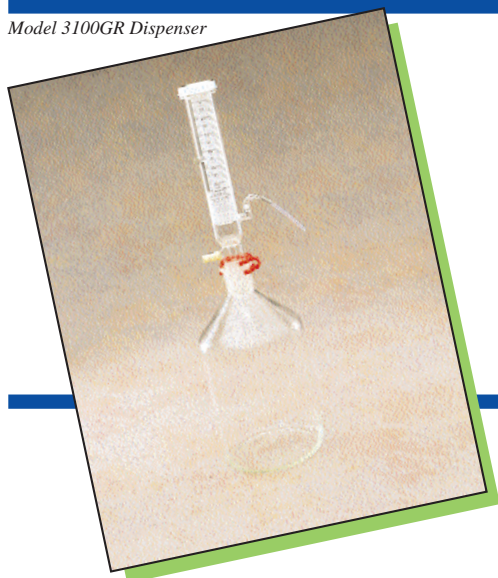

## ALL GLASS DISPENSERS!

- — Ground glass joints
- — Use with an Erlenmeyer flask on your hot plate
- ★ — Glassware carries a 50% lifetime guarantee against breakage

## ORDERING INFORMATION

Model #	Description	Size (ml)	T Joint Size	US List Price
3005GR	Dispenser with 500 ml Reagent Bottle, glass intake tube	5.0 ml	24/40	\$333.00
3010GR	Dispenser with 1000 ml Reagent Bottle, glass intake tube	10.0 ml	29/42	\$356.00
3020GR	Dispenser with 1000 ml Reagent Bottle, glass intake tube	20.0 ml	29/42	\$402.00
3050GR	Dispenser with 2000 ml Reagent Bottle, glass intake tube	50.0 ml	29/42	\$532.00
3100GR	Dispenser with 2000 ml Reagent Bottle, glass intake tube	100.0 ml	29/42	\$571.00
3001/2GV	Dispenser with 100 ml Conical Vessel, glass intake tube	0.5 ml	24/40	\$362.00
3001GV	Dispenser with 100 ml Conical Vessel, glass intake tube	1.0 ml	24/40	\$362.00
3001/2GE	Dispenser with 125 ml Erlenmeyer Flask, glass intake tube	0.5 ml	24/40	\$314.00
3001GE	Dispenser with 125 ml Erlenmeyer Flask, glass intake tube	1.0 ml	24/40	\$314.00
3001/2G	Dispenser-only, Teflon intake tube, use with catalog #5000125 or 1170	0.5 ml	24/40	\$282.00
3001G	Dispenser-only, Teflon intake tube, use with catalog #5000125 or 1170	1.0 ml	24/40	\$282.00
3005G	Dispenser-only, Teflon intake tube, use with catalog #1500500	5.0 ml	24/40	\$306.00
3010G	Dispenser-only, Teflon intake tube, use with catalog #15001000	10.0 ml	29/42	\$320.00
3020G	Dispenser-only, Teflon intake tube, use with catalog #15001000	20.0 ml	29/42	\$366.00
3050G	Dispenser-only, Teflon intake tube, use with catalog #15002000	50.0 ml	29/42	\$449.00
3100G	Dispenser-only, Teflon intake tube, use with catalog #15002000	100.0 ml	29/42	\$489.00

## ACCESSORIES

Catalog #	Description	US List Price
1170	100 ml Conical Vessel	\$78.05
5000125	125 ml Erlenmeyer Flask	\$42.65
1500500	500 ml Corning Bottle	\$36.55
15001000	1000 ml Corning Bottle	\$46.35
15002000	2000 ml Corning Bottle	\$95.40

## PRODUCT DESCRIPTION

- A Standard REPIPET dispenser with standard taper, ground glass joints.
- Full range adjustability, 100-Division Scale.
- 1% accuracy and 0.1% reproducibility at full scale guaranteed.
- Includes both magnifying indicator and metal indicator.
- Autoclavable.

48 Hour or Sooner  
Express Shipping  
Guaranteed!\*




\*see inside front cover for details

# Dispensers/Dilutors

REPIPET II® DISPENSERS

Labindustries

Model 9005 Dispenser

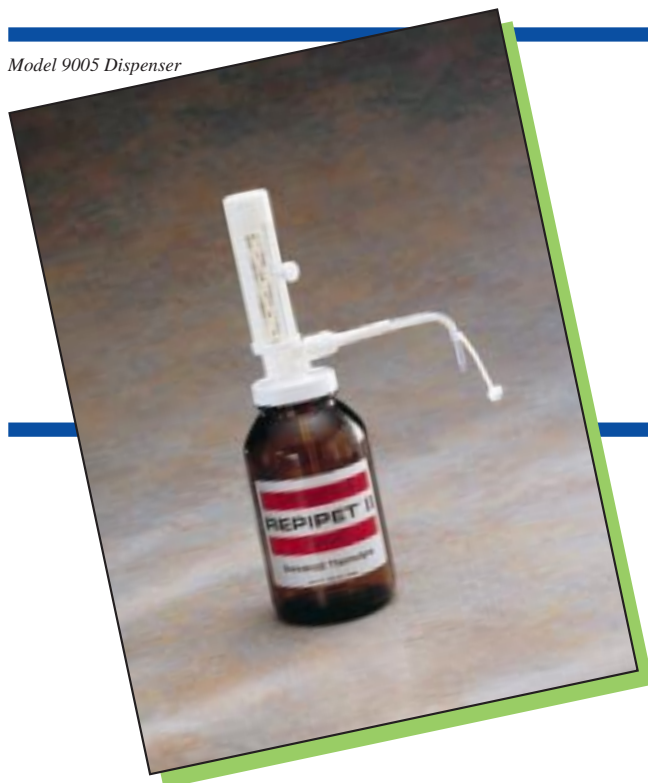

## VIRTUALLY UNBREAKABLE!

- — Magnifying indicator for easy-reading accuracy
- — Integrated 100-division scale
- — Extra thick borosilicate glass barrel and plunger
- ★ — Glassware carries a 50% lifetime guarantee against breakage

## ORDERING INFORMATION

Model #	Size (ml)	Bottle Size	US List Price
Standard REPIPET® II, with amberized glass bottle			
9005	5.0 ml	500 ml	\$132.00
9010	10.0 ml	1250 ml	\$132.00
9020	20.0 ml	250 ml	\$153.00
Universal REPIPET® II, with 3 caps/liners (sizes 28-430, 33-430, 38-430) and 12" (30 cm) Teflon® tubing			
9005U	5.0 ml		\$141.00
9010U	10.0 ml		\$141.00
9020U	20.0 ml		\$161.00

48 Hour or Sooner  
Express Shipping  
Guaranteed!\*



\*see inside front cover for details

# Dispensers/Dilutors

FIXED AND ADJUSTABLE VOLUME

Labindustries

Model 770200 Fixed Volume Dispenser



Includes  
FREE  
Support stand

## OUR LOWEST PRICED DISPENSERS!

- — Fixed volume available in 13 sizes
- — Adjustable models are economical, plastic and unbreakable

## ORDERING INFORMATION

### Fixed Volume Economy

Model #	Size (ml)	Container Size	US List Price
770100	0.1 ml	4 oz. (120 ml)	\$87.00
770200	0.2 ml	4 oz. (120 ml)	\$87.00
770500	0.5 ml	4 oz. (120 ml)	\$87.00
770700	0.7 ml	4 oz. (120 ml)	\$87.00
770900	0.9 ml	8 oz. (240 ml)	\$87.00
7701	1.0 ml	8 oz. (240 ml)	\$87.00
77011/2	1.5 ml	8 oz. (240 ml)	\$87.00
7702	2.0 ml	8 oz. (240 ml)	\$87.00
77021/2	2.5 ml	8 oz. (240 ml)	\$87.00
7703	3.0 ml	8 oz. (240 ml)	\$87.00
7704	4.0 ml	16 oz. (480 ml)	\$87.00
7705	5.0 ml	16 oz. (480 ml)	\$87.00
7710	10.0 ml	16 oz. (480 ml)	\$87.00

### Adjustable Volume Economy

70021	2.0 ml	\$86.00
7005	5.0 ml	\$86.00
7010	10.0 ml	\$86.00

## FIXED VOLUME ECONOMY

- Specially designed recirculating system for priming without loss of reagent.
- Choose any of 13 sizes.
- Guaranteed 2% accuracy and 5% reproducibility.
- Each model includes one amberized glass bottle and one plastic container.
- Adaptable for use with any reagent bottle having a 20-400, 22-400, 24-400, 28-400 or 38-400 screw cap.
- 5 interchangeable cap/adaptor sets for container flexibility.
- Support stand provides added stability to 4 oz. (120 ml) bottle.
- Twelve extra inches (30.5 cm) of Teflon intake tubing.

## ADJUSTABLE VOLUME ECONOMY

- Economical and durable high density polyethylene and polypropylene dispenser.
- Time saving, semi-automatic reloading feature for repetitive lab tests.
- Most economical dispenser on the market today.
- Excellent choice for aqueous solutions and a wide variety of chemicals. (Not suitable for use with oils, hydrocarbons, iodine, petroleum-based chemicals, oxidizing acids, alicyclics or solvents.)

Model 7010 Adjustable Volume Dispenser



Includes  
FREE  
Support stand

48 Hour or Sooner  
Express Shipping  
Guaranteed!\*



\*see inside front cover for details

# Dispensers/Dilutors

## REPIPET® SINGLE STROKE DILUTORS

Labindustries

Model S81001U Dilutor



### COMPARE THE PRICE TO AUTOMATIC DILUTORS!

- — A specialized instrument for diluting samples
- — Dilutes with a single action
- — Adapts to your container
- ★ — Glassware carries a 50% lifetime guarantee against breakage

### ORDERING INFORMATION

Model #	Dilutor	Aspirator	Aspirated Volumes	US List Price
<b>Low Silhouette</b>				
S80101/2L	1.0 ml	0.5 ml	1:1	\$527.00
S81001/2L	10.0 ml	0.75 ml	12:1	\$527.00
S81001L	10.0 ml	4.3 ml	2.3:1	\$527.00
<b>Universal</b>				
S80101/2U	1.0 ml	0.5 ml	1:1	\$508.00
S81001/2U	10.0 ml	0.75 ml	13:1	\$508.00
S81001U	10.0 ml	4.3 ml	2.3:1	\$508.00

### PRODUCT DESCRIPTION

- As fast and precise as expensive automatic dilutors.
- Dilutes with a single action; aspirating and diluting with just one cycle.
- Coupled together by a magnifying indicator coupling; diluent plunger and aspirator plunger raise simultaneously.
- 100-Division Scale.
- Low Silhouette model includes a low-profile, 1000 ml square amberized bottle.
- Universal REPIPET dilutor supplied with three cap and washer sets and 24 inches of Teflon® tubing.
- Both models feature magnifying indicators for easy reading accuracy.
- Pyrex® glass and Teflon® construction make REPIPET dilutors safe for use with all reagents except hydrofluoric solutions. Two-opening, Halar® lids may distort when used with organic solvents.

48 Hour or Sooner  
Express Shipping  
Guaranteed!\*




\*see inside front cover for details

# Dispensers/Dilutors

## AIR DISPLACEMENT PIPETS

Labindustries

Model T705B Pipet

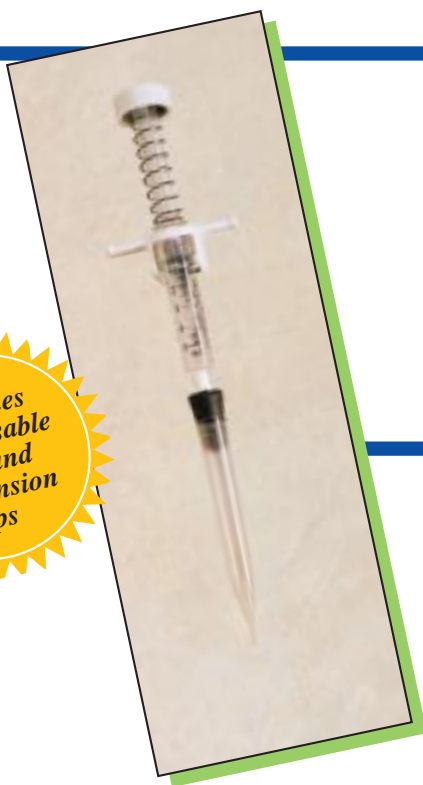

Includes  
3 disposable  
tips and  
3 extension  
tips

### CHEMICALLY RESISTANT!

- — All plastic construction
- — Spring return, syringe type operation

### ORDERING INFORMATION

#### Air Displacement Pipets

Model #	Description	US List Price
T705	Sampling System with stand	\$80.00
T710	10 ml Sampling System with stand	\$80.00
T705B	5 ml Sampling System only	\$64.00
T710B	10 ml Sampling System only	\$64.00

### ACCESSORIES

#### Polypropylene Disposable Tips

Model #	Size ml	Qty/Pkg	US List Price
D7055	5.0	50	\$19.65
D7051	5.0	100	\$39.00
D7105	10.0	50	\$28.60
D7101	10.0	100	\$42.70

### PRODUCT DESCRIPTION

- Sturdy all-plastic construction makes it tough.
- Calibrated in 0.1 ml divisions, this macro transfer pipet is adjustable the full length of the scale.
- Accuracy within  $\pm 2\%$  at full volume. Reproducibility  $\pm 0.5\%$  at full volume.
- Resistant to a wide variety of chemicals. Not to be used with oils, hydrocarbons, toluene, petroleum based chemicals, oxidizing and concentrated acids, aromatic solvents, alicyclics, mercaptol reagents or iodine.
- Supplied with three disposable tips and three extension tips.

48 Hour or Sooner  
Express Shipping  
Guaranteed!\*



\*see inside front cover for details