


# Seropette<sup>®</sup> Digital Pipettors

Light and Comfortable Design  
with Superior Accuracy.

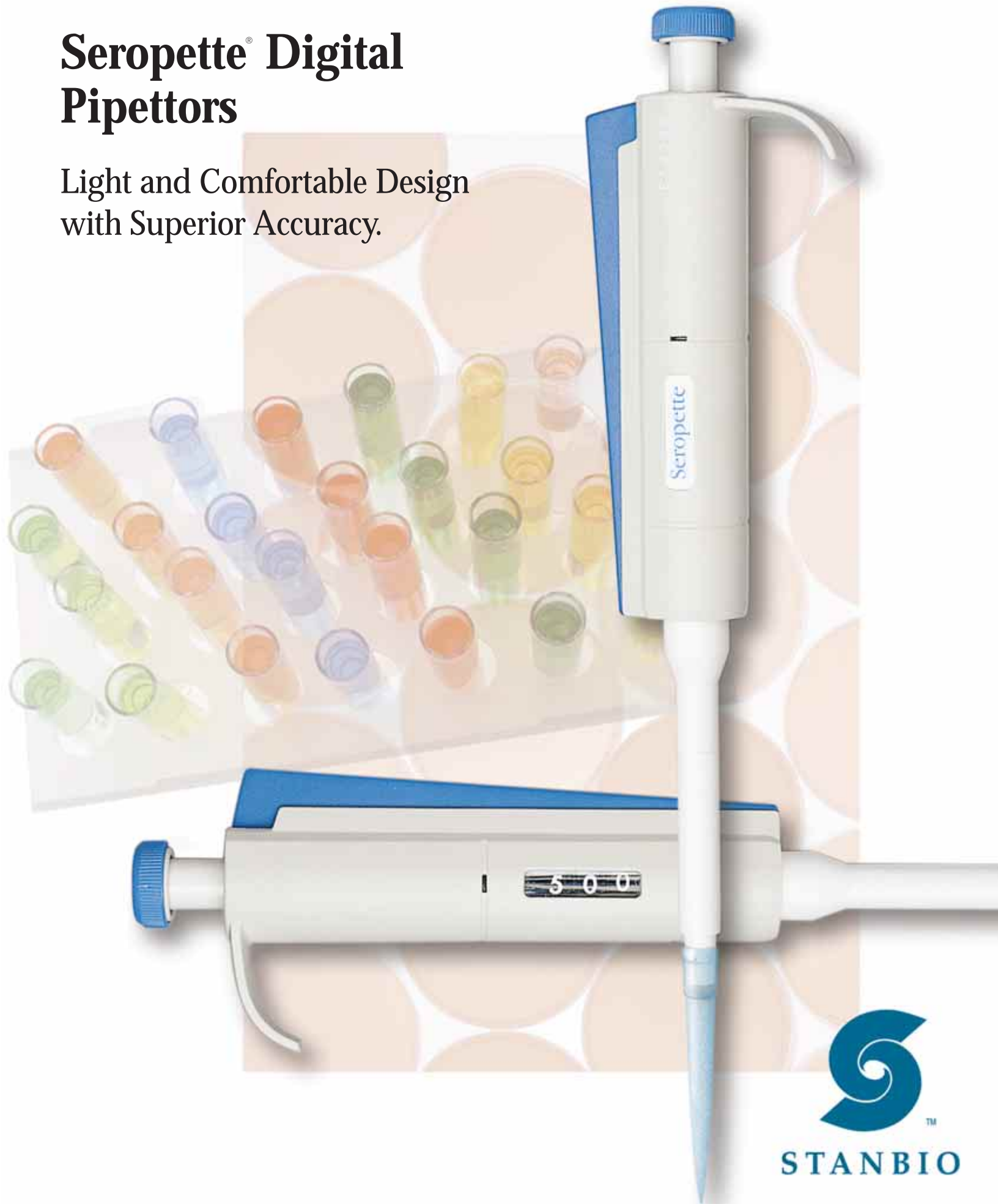

STANBIO

## Seropette® Pipettors – A complete line of variable and fixed-volume models

Light, comfortable design, fixed-volume and variable-volume models with built-in finger rest for minimal hand fatigue. Seropette® pipettors have color-coded control buttons for easy identification of proper pipettor tips. Example: Yellow control button = yellow tip, Blue control button = blue tip.

Built-in tip ejector is separate from

control button to eliminate accidental tip removal. Each pipettor is supplied with an individual calibration certificate to guarantee the highest level of accuracy and precision. All models can be recalibrated in the lab easily with the tool supplied. The tip cone and ejector are fully autoclavable. All models come with a one year warranty and are manufac-

tured under strict ISO 9000 standards. Using Seropette® brand color-coded, disposable pipettor tips made from extra virgin polypropylene assure reliable and precise results. For easy access and proper storage, a durable pipettor stand capable of holding up to 10 Seropette® pipettors is available.

### Seropette® Pipettors



#### Seropette®-Ultra Smooth Laboratory Pipettors

- Comfortable design reduces hand fatigue
- High precision and accuracy - calibration certificate included with every pipettor
- Both fixed-volume and variable-volume available
- Built-in tip ejector for added safety and convenience
- Rugged construction - 1 year warranty
- All pipettors are color - coded for easy selection of Seropette® disposable tips
- Simple to calibrate
- Manufactured following ISO 9000 standards

#### Seropette® Variable-Volume Pipettors

Catalog No.	Volume	Tip Color	Increments	Precision	Accuracy
9000 - 001	5 - 40 µL	Yellow	0.5 µL	2.2 to 0.5 %	± 2.5 to 0.5
9100 - 001	40 - 200 µL	Yellow	1.0 µL	0.8 to 0.4 %	± 1.1 to 0.4
9200 - 001	200 - 1000 µL	Blue	5.0 µL	0.6 to 0.3 %	± 0.9 to 0.3
9300 - 001	1000 - 5000 µL	Natural	5.0 µL	0.6 to 0.3 %	± 0.9 to 0.3

#### Seropette® Fixed-Volume Pipettors

Catalog No.	Volume	Tip Color	Precision	Accuracy
8005 - 001	5 µL	Yellow	2.2%	± 2.5%
8010 - 001	10 µL	Yellow	2.0%	± 1.2%
8020 - 001	20 µL	Yellow	2.0%	± 1.2%
8025 - 001	25 µL	Yellow	1.0%	± 1.0%
8050 - 001	50 µL	Yellow	0.8%	± 1.0%
8100 - 001	100 µL	Yellow	0.6%	± 0.6%
8200 - 001	200 µL	Yellow	0.4%	± 0.4%
8250 - 001	250 µL	Blue	0.6%	± 0.8%
8500 - 001	500 µL	Blue	0.4%	± 0.6%
8001 - 001	1000 µL	Blue	0.3%	± 0.3%

#### Seropette® Disposable Pipettor Tips

Catalog No.	Tip Volume/Color	Tips per Bag	Bags per Case	Total Tips
9500 - 001	5 - 200 µL / Yellow	500 tips	20 bags	10,000
9600 - 001	200 - 1000 µL / Blue	250 tips	40 bags	10,000
9700 - 001	1000 - 5000 µL / Natural	500 tips	1 bag	500

#### Seropette® Stand

Catalog No. 9900-001  
Holds up to 10 Seropette® Pipettors.

**DISTRIBUTED BY:**

*To order, or for more information, call U.S. Toll-Free 1-800-531-5535 or (210) 222-2108*  
Stanbio Laboratory, Inc. • 2930 East Houston Street • San Antonio, Texas 78202 • Fax: (210) 227-6367 or visit our website at: [www.stanbio.com](http://www.stanbio.com)



**STANBIO**

*Advancing the Science of Clinical Diagnostics*

General & Special Chemistry / Pregnancy / Serology / Microbiology / Hematology / Occult Blood / Urinalysis / Blood Grouping / Clinical Instrumentation / Pipettors