

NEW

Easyplate™ Petri Dish Racks

Perfect for Use in an Incubator.

Increase workspace and efficiently transport large quantities of Petri dishes with the Easyplate™ Petri Dish Rack. Silicone holding strips on the rack's four columns grip each Petri dish and lid to enable the user to easily insert and remove groups of dishes without disturbing upper or lower dishes. Labels are provided on the sides of the rack for critical sample identification. Available in 4 colors and two sizes: the 9" size holds forty-eight 100mm Petri Dishes and the 14" size holds eighty-four 100mm Petri Dishes. The Easyplate™ rack is manufactured from lightweight, chemically resistant polycarbonate and is autoclavable at 121° C (250°F).

CATALOG NO.	COLOR	SIZE HEIGHT	PRICE/EACH
F18982-0014	Natural	369mm (14 ⁹ / ₁₆ ")	\$99.50
F18982-0314	Blue	369mm (14 ⁹ / ₁₆ ")	99.50
F18982-0009	Natural	240mm (9 ¹¹ / ₁₆ ")	92.50
F18982-0109	Red	240mm (9 ¹¹ / ₁₆ ")	92.50
F18982-0209	Yellow	240mm (9 ¹¹ / ₁₆ ")	92.50
F18982-0309	Blue	240mm (9 ¹¹ / ₁₆ ")	92.50
F18982-0409	Natural, Red, Yellow, Blue	240mm (9 ¹¹ / ₁₆ ")	355.00/Set 4

Easyplate™ Turntable

Using the Easyplate™ Turntable the researcher or technician can rotate the Easyplate™ Petri Dish Rack 360° for fast, easy accessibility. The turntable revolves smoothly on ball bearings, allowing the user to quickly remove or replace a Petri dish on any side of the Easyplate™ Petri Dish Rack. The 6-inch post in the center, made of high-density polyethylene, orients and stabilizes the rack. The turntable is white styrene and is easily cleaned with alcohol-based disinfectant.

CATALOG NO.	PRICE/EACH
F18982-1000	\$31.00

NEW

100mm Petri Dish Rack

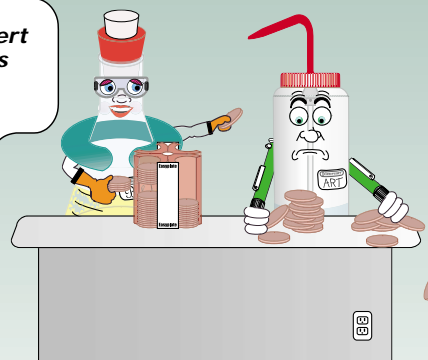
Maximize Incubator Space

Use this acrylic rack to hold up to fifty-four 100mm Petri dishes in organized stacks of 9 dishes. Dishes can be removed with one hand and are easily viewed through the clear base and divider. The divider has a built in finger grip handle for transporting the rack when loaded with Petri dishes. Rubber feet help prevent sliding on the bench

top. Racks are shipped flat and can be assembled in minutes. Dimensions: 34L X 20W X 23cm H (13¹/₂ L X 9¹/₄ H X 8" W)

CATALOG NO.	PRICE/EACH
F18985-0100	70.00

Remove and insert multiple dishes at once.

**NEW**

Finger slot for easy removal of dishes

Molded 100mm Petri Dish Rack

Safely Transports Forty-Two Dishes

This rack is useful for safe transport and incubation of 100mm Petri dishes and helps to streamline benchtop sample processing. Each rack holds up to forty-two Petri dishes with lids, configured in six columns of seven plates. The base area is large enough to be stable but small enough to preserve

bench top space. Up to four racks can be stacked together, and all interlock to prevent an accidental knock over. Access slots under the bottom Petri dishes allows a stack of dishes to be safely removed from the rack using a secure grip that places a finger under the bottom dish and thumb on top of the lid of the top dish. A large center divider provides space for labeling and incorporates a handle that does not interfere with the samples. To maintain sterility and ensure long term performance, all parts are molded with white polypropylene and are steam autoclavable at 121°C (250°F). Racks are shipped flat and can be easily assembled. Overall dimensions are 33L x 21W x 17.8H cm (13 L x 8¹/₄ W x 7" H)

CATALOG NO.	PRICE/EACH
F18991-0000	\$39.00

Petri Dish Dispenser Rack

Adjusts to Hold a Variety of Dish Sizes

This adjustable dispenser neatly and safely holds 15 Petri dishes up to 102mm (4") diameter. Made of rugged ABS plastic, it slides apart at the base for easy refilling. The rack is 89mm (3¹/₂") wide, extends from 127mm (5") to 210mm (8³/₄") long and is 210mm (8³/₄") high.

CATALOG NO.	PRICE/EACH
F18982-0002	\$34.50

Poxygrid® Petri Dish Mini Rack

Organize and Safely Hold Six Dishes

This Petri Dish Mini Rack helps keep benchtops neat and organized. Each welded steel wire rack conveniently holds six 100mm Petri dishes with lids. The rack is ideal for lab transport and incubation of dishes. Durable blue epoxy coated racks are steam autoclavable at 121°C (250°F). Overall dimensions are 13 x 8.5 x 12H cm (5¹/₄ x 3¹/₄ x 4³/₄" H).



CATALOG NO.	PRICE/EACH
F18979-0006	\$11.95